


**NovaPED** calcaneus  
by **schein**

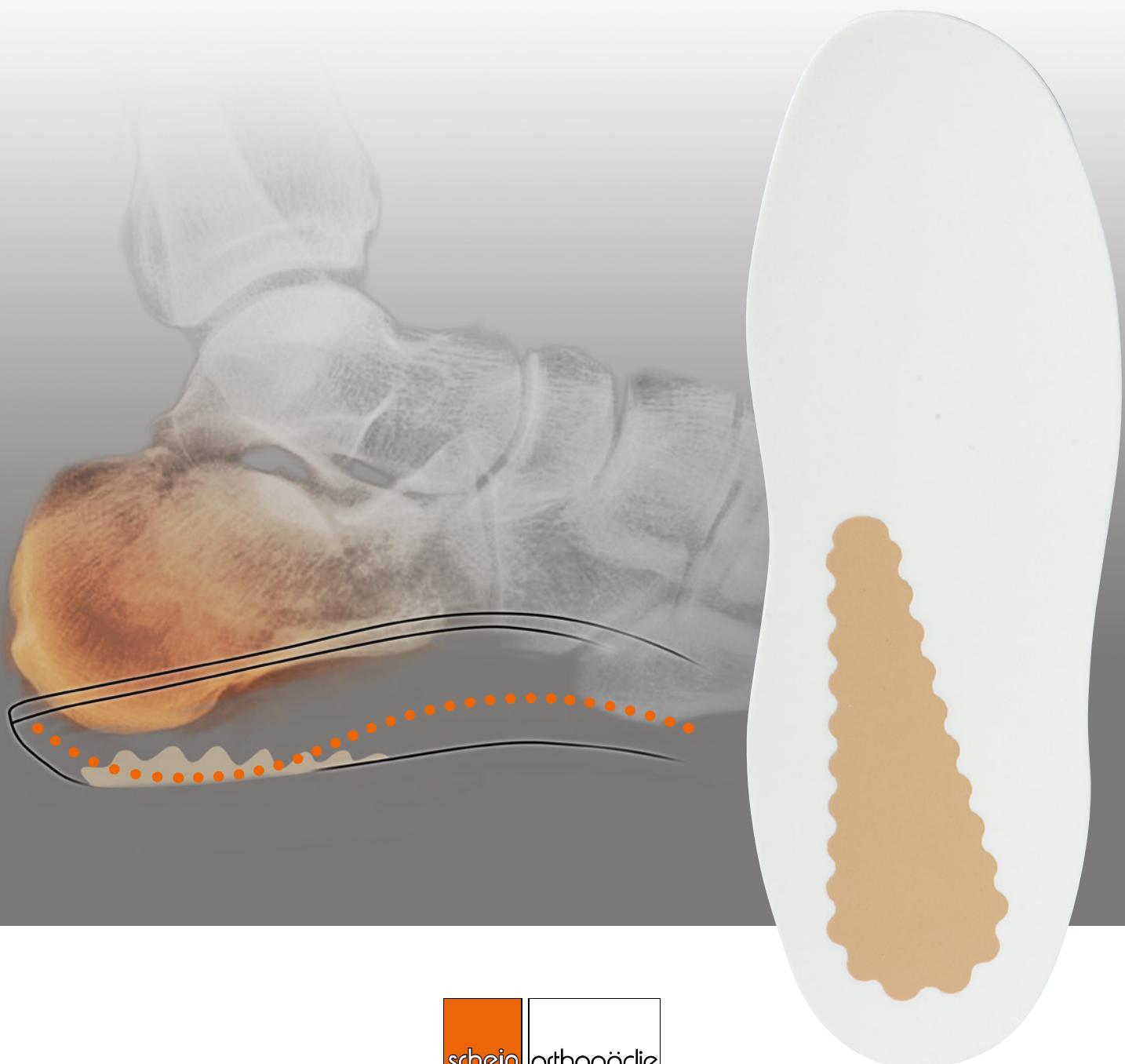

**schein** orthopädie service  
Germany since 1879

WeXhelp

INCREASING QUALITY OF LIFE  
AND WELLBEING OF PEOPLE!

## Novaped S90 calcaneus

### NOVAPED S90 CALCANEUS BUSINESS

Art. no.	Description	Size
989818000	Novaped S90 calcaneus business	35–50
989818174	Novaped S90 calcaneus business, loose cover	35–50



### NOVAPED S90 CALCANEUS SNEAKER

Art. no.	Description	Size
989820000	Novaped S90 calcaneus Sneaker	35–50
989820174	Novaped S90 calcaneus Sneaker, loose cover	35–50



## Novaped soft calcaneus

### NOVAPED SOFT CALCANEUS RESORB SOFT

Art. no.	Description	Size
979828000	Novaped soft calcaneus resorb soft	35–50



### NOVAPED SOFT CACANEUS ABSORB MEDIUM

Art. no.	Description	Size
979825000	Novaped soft calcaneus absorb medium	35–50



# NovaPED calcaneus

by schein

## NOVAPED CALCANEUS


The **NovaPED calcaneus** series has been specially designed for use in case of overload symptoms or pain in the area of the calcaneus/plantar tendon attachment.

The **NovaPED S90 calcaneus** orthopaedic foot supports are preformed from a thermoplastic material and, in addition to an anatomically adapted recess in the middle of the heel, also have a lowering of the plantar tendon structure.

The use of glass-fibre-reinforced components increases the stability in the longitudinal arch, while PU padding materials in the heel have a shock-absorbing and pressure-relieving effect. Another special feature is the multi-zone technology in the cover combination. The cushioning top layer with soft-bedding heel inlay supports and relieves the foot structure in a targeted way.

The multi-zone technology is also used in the **NovaPED soft calcaneus** orthopaedic foot supports. The base made of Multiforte provides the form and is functionally equipped through the use of pressure-relieving or absorbent PU materials.

The blanks of the NovaPED calcaneus series can be thermally formed/shaped/milled and are intended for further processing by a healthcare technician.

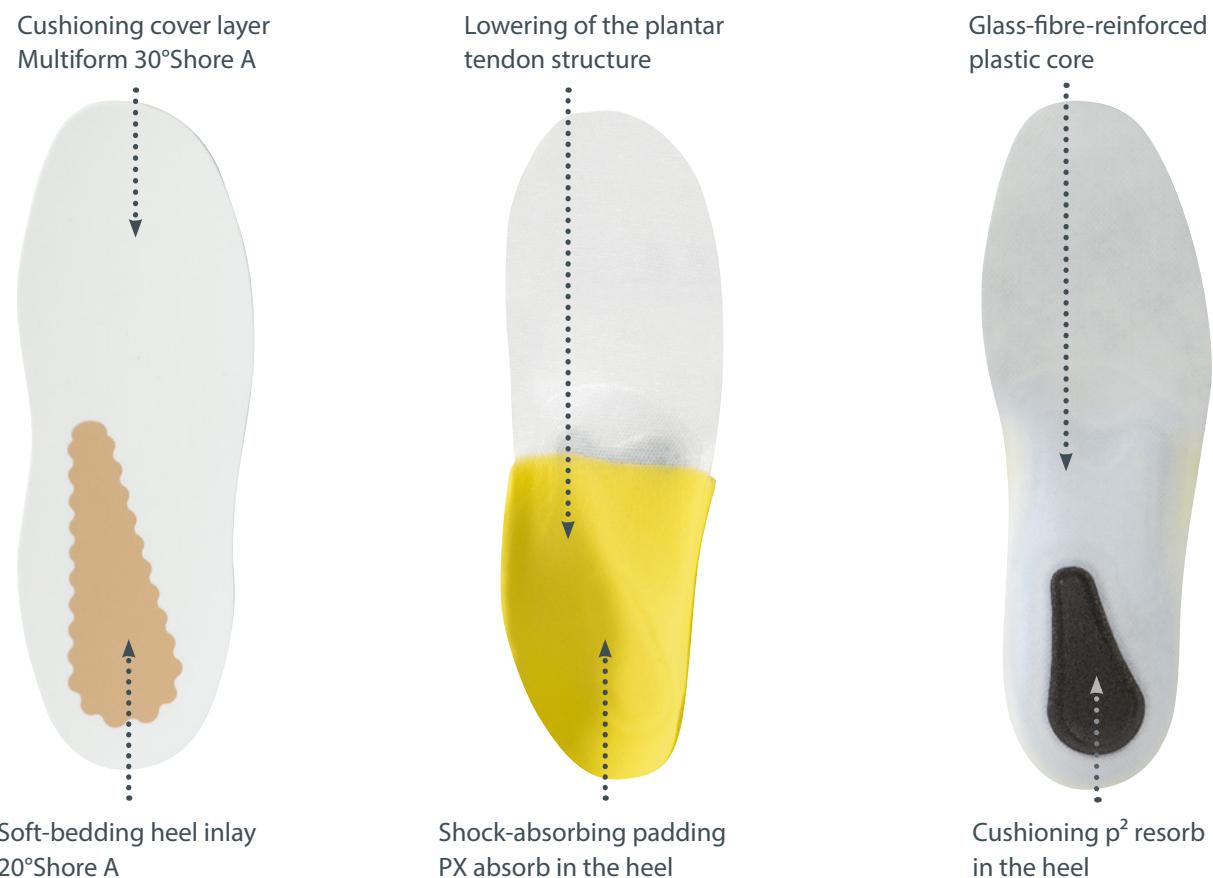

# Novaped S90 calcaneus

## NOVAPED S90 CALCANEUS BUSINESS

- Supporting insoles produced using the film technique
- Glass-fibre-reinforced plastic core
- Lowering of the plantar tendon structure
- Heel padding made of PU cushioning elements (p<sup>2</sup> resorb/PX absorb)
- Anatomical rubber padding
- OnDry cover, dark grey, embossed
- Can be thermally moulded

Especially thin  
in the forefoot area



Art. no.	Description	Size
989818000	Novaped S90 calcaneus business	35-50
989818174	Novaped S90 calcaneus business, loose cover	35-50

## NOVAPED S90 CALCANEUS SNEAKER

- Supporting insoles produced using the film technique
- Glass-fibre-reinforced plastic core
- Lowering of the plantar tendon structure
- Heel padding made of PU cushioning elements (p<sup>2</sup> resorb/PX absorb)
- Anatomical rubber padding
- Long-soled Multiform padding in a multi-zone technique, 3 mm
- OnDry cover, dark grey, embossed
- Can be thermally moulded



Art. no.	Description	Size
989820000	Novaped S90 calcaneus Sneaker	35-50
989820174	Novaped S90 calcaneus Sneaker, loose cover	35-50

# Novaped soft calcaneus

## NOVAPED SOFT CALCANEUS RESORB SOFT

- Bedding orthopaedic foot supports produced using soft cushioning technique ( $p^2$  resorb)
- Base in Multiforte multi-zone technique (Multiforte/ $p^2$  resorb)
- Especially soft heel impact zone
- Reinforcing intermediate layer made of THK (windowed)
- Shaped soft pad
- Long-soled Multiform padding in a multi-zone technique, 4 mm
- OnDry cover, dark grey, embossed
- Can be thermally shaped



Art. no.	Description	Size
979828000	Novaped soft calcaneus resorb soft	35–50

## NOVAPED SOFT CALCANEUS ABSORB MEDIUM

- Bedding orthopaedic foot supports produced using soft cushioning technique ( $p^2$  absorb)
- Base in Multiforte multi-zone technique (Multiforte/ $p^2$  absorb)
- Energy-absorbing heel padding
- Reinforcing intermediate layer made of THK (windowed)
- Shaped soft pad
- Long-soled Multiform padding in a multi-zone technique, 4 mm
- OnDry cover, dark grey, embossed
- Can be thermally shaped



Art. no.	Description	Size
979825000	Novaped soft calcaneus absorb medium	35–50



099993217

Schein Orthopädie Service KG  
Hildegardstr. 5  
42897 Remscheid  
Germany  
Tel. +49 2191 910-200  
Fax +49 2191 910-100  
[service@schein.de](mailto:service@schein.de)  
[www.schein.de](http://www.schein.de)



orthopädie  
service

— Germany since 1879 —

We help