



NovaPED

by **schein**

soft



Einlagenrohlinge aus Weichpolstermaterial haben in den vergangenen Jahren eine dominierende Rolle übernommen. Dank der zahlreichen Variationsmöglichkeiten hinsichtlich Härte, Flexibilität, Thermoplastizität, Rückstellfähigkeit und nicht zuletzt auch der Farbe kann praktisch für jeden Fuß der passende Rohling erstellt werden. Voraussetzung für eine optimale Versorgung ist natürlich auch die Kenntnis bei den Leistungserbringern über die verschiedenen Qualitäten der Weichpolstermaterialien. Dass heute Kunststoffe bedarfsgerecht hergestellt werden können, ist ein großer Fortschritt der chemischen Industrie, der letztlich den Patienten mit Fußbeschwerden zugute kommen sollte, der andererseits aber auch die stete Information voraussetzt. Seit Jahrzehnten wird bei Schein an immer wieder neuen Materialien experimentiert und Anregungen aus der Praxis auf Einsatzmöglichkeit für die industrielle Fertigung überprüft. So konnte manche Neuentwicklung auf den Weg gebracht werden. Erfreulich für die Anwender ist die Erkenntnis, dass es nur wenige Grundstoffe bei den Weichpolstermaterialien gibt, die sich jedoch durch die genannten Qualitäten unterscheiden. Deshalb soll eine kurze Charakteristik zur Gebrauchsinformation dienen:

Der am meisten eingesetzte Kunststoffschaum besteht aus **Polyethylen** (auch Polyethen oder **PE** genannt). Er ist reiner Kohlenwasserstoff, der sich aus vielen Ethylenmolekülen ($\text{CH}_2 = \text{CH}_2$) durch Polymerisation zusammenfügt und ohne Rückstände entsorgt werden kann. Er ist in folgenden von uns eingesetzten Materialien zu finden:

Multiform, Multiforte, PE-Schaum grau.

Vorteile von PE:

thermisch verformbar, sehr viele Shore-Grade und Farben verfügbar.

Ein Material, was in den letzten Jahren explosionsartig den Markt bereichert hat, ist **EthylVinylAcetat**, besser als **EVA** bekannt. Es handelt sich dabei um eine Gruppe von Copolymeren, die aus Ethylen, Vinyl und Acetat in unterschiedlichen Molekülmengen zusammengefügt werden. EVA findet sich in folgenden unserer Weichpolstermaterialien: **Pryx, Thermofix, Pryx antistatisch und natürlich EVA.**

Vorteile von EVA:

thermisch verformbar, sehr viele Shore-Grade (18°- 70°), sehr gute Rückstellfähigkeit, problemlose Be- und Verarbeitung, auch antistatisch in Sicherheitsschuhen verwendbar.

Als dritter Grundstoff soll das **Polyurethan (PU oder PUR)** beschrieben werden. Es ist ein Polyaddukt aus Grundmolekülen von Polyisocyanaten und Polyalkoholen. Durch die innere Vernetzung der Moleküle werden die weichen Polyurethane den Elastomeren zugeordnet, d.h., dass sie zwar elastisch aber nicht thermoplastisch sind. Bei unseren Einlagenrohlingen befindet sich Polyurethan in folgenden Weichpolstermaterialien: **ARUfoam, PU hautfarben, Multicom (PE + PU) und Plastocom (PE + PU).**

Vorteile von PUR:

hohe Rückstellfähigkeit, offene Poren, Verschleißfestigkeit, bei der Kombination mit PE auch thermoplastisch.

During recent years, soft-padded orthopaedic foot supports have adopted a commanding role. Thanks to the numerous variations regarding hardness, flexibility, thermoplasticity, resilience and last but not least also regarding colour, it is possible today to make an appropriate foot orthosis for every foot. Of course, a prerequisite for optimum care is that the care providers have appropriate knowledge about the various qualities of the soft-padded materials. The fact that plastic materials can be produced today which meet our demands is a huge progress of the chemical industry, which in the end should be of benefit to the patient with foot disorders, however, on the other hand, this also presupposes continuous information. We at Schein have been experimenting with new materials for decades and have checked new ideas from clinical practice for possible applications for industrial manufacture. This way, many a new developments could be launched on the market. Enjoyable to the users is the awareness that there are only few base materials among the soft-padded materials which, however, are distinguished by the qualities mentioned before. Therefore, a short characterization should provide some information on use:

*The most frequently used plastic material consists of **polyethylene** (also called polyethene or **PE**). This is pure hydrocarbon which merges by polymerisation from many ethylene molecules ($\text{CH}_2 = \text{CH}_2$) and which can be disposed of without residues. It is contained in the following materials that we use:*

Multiform, Multiforte, PE-Foam grey.

PE benefits:

thermoformable, available in a great many Shore grades and colours.

*Another material, which during recent years has conquered the market with explosive growth, is **EthylVinylAcetate**, better known as **EVA**. This material can be described as a group of copolymers that are joined together from ethylene, vinyl and acetate in different molecular amounts. EVA is used in the following of our soft-padded materials:*

Pryx, Thermofix, Pryx anti-static and, of course, EVA.

EVA benefits:

thermoformable, a great many Shore grades (18°- 70°), extremely good resilience, trouble-free treatment and processing, and it can also be used as an antistatic material in safety shoes.

*A third base material to be described here is the **polyurethane (PU or PUR)**. This is a polyadduct of base molecules from polyisocyanates and polyalcohols. Due to the inner cross-linking of the molecules the soft polyurethanes are allocated to the elastomers, i.e. that they are elastic indeed, but not thermoplastic. We use polyurethane for our orthopaedic foot supports in the following soft-padded materials:*

ARUfoam, skin-coloured PU, Multicom (PE + PU) and Plastocom (PE + PU).

PUR benefits:

high resilience, open pores, wear resistance, when combined with PE also thermoplastic.

NovaPED soft

Seite/page

Einschicht / *single layer*

04 - 09

Mehrschicht / *multi layer*

10 - 19

Material <i>Material</i>	Chemischer Grundstoff <i>chemical base material</i>	ca. Shore A <i>approx. density</i>	Materialtemperatur für Verformung <i>working temperature</i>	Empfohlener Kleber <i>recommended glue</i>	Vorschleiff <i>pre-finishing</i>	Fertigschleiff <i>finishing</i>
PE-Schaum, grau <i>PE-Foam, grey</i>	PE	20	max. 90 - 100° C	Renia Ortec Spezial	Korn 40, Schleifband <i>grain 40, sanding belt</i>	Korn 80, Schleifband <i>grain 80, sanding belt</i>
Multiform <i>Multiform</i>	PE / EVA	30	max. 90 - 100° C	Renia Ortec Spezial	Korn 40, Schleifband <i>grain 40, sanding belt</i>	Korn 80, Schleifband <i>grain 80, sanding belt</i>
ARUform 20 <i>ARUform 20</i>	PE	36	max. 90 - 100° C	Renia Ortec Spezial	Korn 40, Schleifband <i>grain 40, sanding belt</i>	Korn 80, Schleifband <i>grain 80, sanding belt</i>
Multiforte <i>Multiforte</i>	PE / EVA	40	max. 90 - 100° C	Renia Ortec Spezial	Korn 40, Schleifband <i>grain 40, sanding belt</i>	Korn 80, Schleifband <i>grain 80, sanding belt</i>
Multiform, fest <i>Multiform, hard</i>	PE / EVA	50	max. 90 - 100° C	Renia Ortec Spezial	Korn 40, Schleifband <i>grain 40, sanding belt</i>	Korn 80, Schleifband <i>grain 80, sanding belt</i>
Thermaforte <i>Thermaforte</i>	EVA	70	max. 90 - 100° C	Renia Ortec Spezial	Korn 40, Schleifband <i>grain 40, sanding belt</i>	Korn 80, Schleifband <i>grain 80, sanding belt</i>
ARUped 18 <i>ARUped 18</i>	EVA	18	max. 90 - 100° C	Renia Ortec Spezial	Korn 40, Schleifband <i>grain 40, sanding belt</i>	Korn 80, Schleifband <i>grain 80, sanding belt</i>
ARUped22 <i>ARUped22</i>	EVA	22	max. 90 - 100° C	Renia Ortec Spezial	Korn 40, Schleifband <i>grain 40, sanding belt</i>	Korn 80, Schleifband <i>grain 80, sanding belt</i>
ARUped25 <i>ARUped25</i>	EVA	25	max. 90 - 100° C	Renia Ortec Spezial	Korn 40, Schleifband <i>grain 40, sanding belt</i>	Korn 80, Schleifband <i>grain 80, sanding belt</i>
Pryx 25 <i>Pryx 25</i>	EVA	25	max. 90 - 100° C	Renia Ortec Spezial	Korn 40, Schleifband <i>grain 40, sanding belt</i>	Korn 80, Schleifband <i>grain 80, sanding belt</i>
Thermofix <i>Thermofix</i>	EVA	45	max. 90 - 100° C	Renia Ortec Spezial	Korn 40, Schleifband <i>grain 40, sanding belt</i>	Korn 80, Schleifband <i>grain 80, sanding belt</i>
Pryx, antistatisch <i>Pryx, anti-static</i>	EVA	55	max. 90 - 100° C	Renia Renol Asti Spezial	Korn 40, Schleifband <i>grain 40, sanding belt</i>	Korn 80, Schleifband <i>grain 80, sanding belt</i>
Multikork <i>Multikork</i>	EVA	60	max. 90 - 100° C	Renia Ortec Spezial	Korn 40, Schleifband <i>grain 40, sanding belt</i>	Korn 80, Schleifband <i>grain 80, sanding belt</i>
p ² <i>p²</i>	PU	18		Renia Ortec Spezial	Korn 40, Schleifband <i>grain 40, sanding belt</i>	Korn 80, Schleifband <i>grain 80, sanding belt</i>
x ² <i>x²</i>	PU	15		Renia Ortec Spezial	Korn 40, Schleifband <i>grain 40, sanding belt</i>	Korn 80, Schleifband <i>grain 80, sanding belt</i>
ARUfoam <i>ARUfoam</i>	PU	14		Renia Ortec Spezial	Korn 40, Schleifband <i>grain 40, sanding belt</i>	Korn 80, Schleifband <i>grain 80, sanding belt</i>
Multicom <i>Multicom</i>	PE / EVA / PU	30/18		Renia Ortec Spezial	Korn 40, Schleifband <i>grain 40, sanding belt</i>	Korn 80, Schleifband <i>grain 80, sanding belt</i>
Plastocom <i>Plastocom</i>	PE/PU	20/18		Renia Ortec Spezial	Korn 40, Schleifband <i>grain 40, sanding belt</i>	Korn 80, Schleifband <i>grain 80, sanding belt</i>

Einschicht
single layer


971303-000



971310-000



971303-006



971303-133

Artikel-Nr. Code-No.	Unterbau/Kern Base	Decke Top Layer	Bezug Top Cover	Ausführung Design	Größen Sizes
971303- 000		Multiform, weiss, Shore ca. 30°, 15 mm <i>Multiform, white, Shore approx. 30°, 15 mm</i>		lang, mittel, Pelottendruck <i>full length, medium width, natural meta pad</i>	35 - 48
971310- 000		Multiform, weiss, Shore ca. 30°, 15 mm <i>Multiform, white, Shore approx. 30°, 15 mm</i>		lang, breit, Pelottendruck <i>full length, wide, natural meta pad</i>	19/20 - 50
971303- 006		Multiform, weiss, Shore ca. 30°, 15 mm <i>Multiform, white, Shore approx. 30°, 15 mm</i>		lang, mittel, Pelottendruck, G-schliff <i>full length, medium width, natural meta pad, grinded bottom</i>	35 - 48
971310- 006		Multiform, weiss, Shore ca. 30°, 15 mm <i>Multiform, white, Shore approx. 30°, 15 mm</i>		lang, breit, Pelottendruck, G-schliff <i>full length, wide, natural meta pad, grinded bottom</i>	19/20 - 50
971303- 133		Multiform, weiss, Shore ca. 30°, 15 mm <i>Multiform, white, Shore approx. 30°, 15 mm</i>	Rindleder, perforiert, mild <i>neat's leather, perforated, soft</i>	lang, mittel, Pelottendruck, G-schliff <i>full length, medium width, natural meta pad, grinded bottom</i>	35 - 48
971310- 133		Multiform, weiss, Shore ca. 30°, 15 mm <i>Multiform, white, Shore approx. 30°, 15 mm</i>	Rindleder, perforiert, mild <i>neat's leather, perforated, soft</i>	lang, breit, Pelottendruck <i>full length, wide, natural meta pad</i>	19/20 - 50
971303- 300		Multiform, weiss, Shore ca. 30°, 15 mm <i>Multiform, white, Shore approx. 30°, 15 mm</i>		lang, mittel, ohne Pelottendruck <i>full length, medium width, without meta pad</i>	35 - 48
971310- 300		Multiform, weiss, Shore ca. 30°, 15 mm <i>Multiform, white, Shore approx. 30°, 15 mm</i>		lang, breit, ohne Pelottendruck <i>full length, wide, without meta pad</i>	19/20 - 50

Einschicht
single layer

Artikel-Nr. Code-No.	Unterbau/Kern Base	Decke Top Layer	Bezug Top Cover	Ausführung Design	Größen Sizes
971310-467		Multiform, weiss, Shore ca. 30°, 15 mm	Alcantara, beige	lang, breit, Pelottendruck	19/20 - 50
		Multiform, white, Shore approx. 30°, 15 mm	Alcantara, beige	full length, wide, natural meta pad	
972030-000		Multiform, weiss, Shore ca. 30°, 15 mm		kurz, mittel, Pelottendruck	35 - 48
		Multiform, white, Shore approx. 30°, 15 mm		short, medium width, natural meta pad	
972100-000		Multiform, weiss, Shore ca. 30°, 15 mm		kurz, breit, Pelottendruck	35 - 48
		Multiform, white, Shore approx. 30°, 15 mm		short, wide, natural meta pad	
972030-300		Multiform, weiss, Shore ca. 30°, 15 mm		kurz, mittel, ohne Pelottendruck	35 - 48
		Multiform, white, Shore approx. 30°, 15 mm		short, medium width, without meta pad	
972100-300		Multiform, weiss, Shore ca. 30°, 15 mm		kurz, breit, ohne Pelottendruck	35 - 48
		Multiform, white, Shore approx. 30°, 15 mm		short, wide, without meta pad	



971310-467



972030-000




972100-000

Einschicht
single layer


470770-000



470772-000

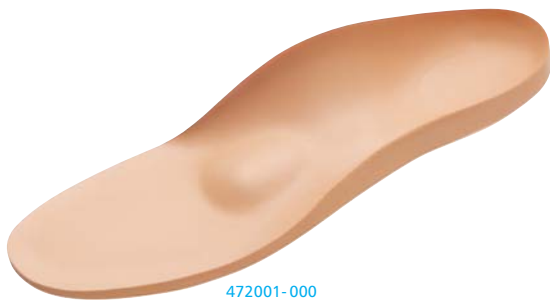

472001-000



472002-000



971320-000



977851-000

Artikel-Nr. Code-No.	Unterbau/Kern Base	Decke Top Layer	Bezug Top Cover	Ausführung Design	Größen Sizes
470770-000		Multiform, weiss, Shore ca. 30°, 15 mm <i>Multiform, white, Shore approx. 30°, 15 mm</i>		lang, breit, Pelottendruck <i>full length, wide, natural meta pad</i>	19/20 - 48
470772-000		Multiform, weiss, Shore ca. 30°, 15 mm <i>Multiform, white, Shore approx. 30°, 15 mm</i>		lang, mittel, Pelottendruck, Unterseite flach <i>full length, medium width, natural meta pad, flat bottom</i>	19/20 - 48
472001-000		Multiform, hautfarben, Shore ca. 30°, 15 mm <i>Multiform, skin-colour, Shore approx. 30°, 15 mm</i>		lang, breit, Pelottendruck <i>full length, wide, natural meta pad</i>	19/20 - 48
472002-000		Multiform, blau, Shore ca. 30°, 15 mm <i>Multiform, blue, Shore approx. 30°, 15 mm</i>		lang, mittel, Pelottendruck <i>full length, medium width, natural meta pad</i>	19/20 - 48
971320-000		Multiform, weiss, Shore ca. 30°, 15 mm <i>Multiform, white, Shore approx. 30°, 15 mm</i>		lang, extrabreit, Pelottendruck <i>full length, extra-wide, natural meta pad</i>	35 - 48
977851-000		Multiform, weiss, Shore ca. 30°, 15 mm <i>Multiform, white, Shore approx. 30°, 15 mm</i>		Naturform, lang, Pelottendruck <i>natural shape, full length, natural meta pad</i>	35 - 48
470770-133		Multiform, weiss, Shore ca. 30°, 15 mm <i>Multiform, white, Shore approx. 30°, 15 mm</i>	Rindleder, perforiert, mild <i>neat's leather, perforated, soft</i>	lang, breit, Pelottendruck <i>full length, wide, natural meta pad</i>	19/20 - 48
971320-133		Multiform, weiss, Shore ca. 30°, 15 mm <i>Multiform, white, Shore approx. 30°, 15 mm</i>	Rindleder, perforiert, mild <i>neat's leather, perforated, soft</i>	lang, extrabreit, Pelottendruck <i>full length, extra-wide, natural meta pad</i>	35 - 48

Einschicht
single layer

Artikel-Nr. Code-No.	Unterbau/Kern Base	Decke Top Layer	Bezug Top Cover	Ausführung Design	Größen Sizes
971320-300		Multiform, weiss, Shore ca. 30°, 15 mm <i>Multiform, white, Shore approx. 30°, 15 mm</i>		lang, extrabreit, ohne Pelottendruck <i>full length, extra-wide, without meta pad</i>	35 - 48
977851-300		Multiform, weiss, Shore ca. 30°, 15 mm <i>Multiform, white, Shore approx. 30°, 15 mm</i>		Naturform, lang, ohne Pelottendruck <i>natural shape, full length, without meta pad</i>	35 - 48



Einschicht
single layer


977013-000



971306-000



971578-000

Artikel-Nr. Code-No.	Unterbau/Kern Base	Decke Top Layer	Bezug Top Cover	Ausführung Design	Größen Sizes
977013-000		Multiform, weiss, Shore ca. 30°, 15 mm <i>Multiform, white, Shore approx. 30°, 15 mm</i>		Viertelschale, lang, Pelottendruck <i>quarter cup, full length, natural meta pad</i>	19/20 - 48
971306-000		Multiform, weiss, Shore ca. 30°, 15 mm <i>Multiform, white, Shore approx. 30°, 15 mm</i>		Schale, lang, Pelottendruck <i>cup, full length, natural meta pad</i>	35 - 48
971578-000		Multiform, weiss, Shore ca. 30°, 15 mm <i>Multiform, white, Shore approx. 30°, 15 mm</i>		Schale, Naturform, lang, Pelottendruck <i>cup, natural shape, full length, natural meta pad</i>	35 - 48
971306-006		Multiform, weiss, Shore ca. 30°, 15 mm <i>Multiform, white, Shore approx. 30°, 15 mm</i>		Schale, lang, Pelottendruck, G-schliff <i>cup, full length, natural meta pad, grinded bottom</i>	35 - 48
977013-006		Multiform, weiss, Shore ca. 30°, 15 mm <i>Multiform, white, Shore approx. 30°, 15 mm</i>		Viertelschale, lang, Pelottendruck, G-schliff <i>quarter cup, full length, natural meta pad, grinded bottom</i>	19/20 - 48
971306-300		Multiform, weiss, Shore ca. 30°, 15 mm <i>Multiform, white, Shore approx. 30°, 15 mm</i>		Schale, lang, ohne Pelottendruck <i>cup, full length, without meta pad</i>	35 - 48
977013-300		Multiform, weiss, Shore ca. 30°, 15 mm <i>Multiform, white, Shore approx. 30°, 15 mm</i>		Viertelschale, lang, ohne Pelottendruck <i>quarter cup, full length, without meta pad</i>	19/20 - 48

Einschicht
single layer

Artikel-Nr. Code-No.	Unterbau/Kern Base	Decke Top Layer	Bezug Top Cover	Ausführung Design	Größen Sizes
470870-000		Multiforte, weiss, Shore ca. 40°, 10 mm <i>Multiforte, white, Shore approx. 40°, 10 mm</i>		kurz, mittel, Pelottendruck <i>short, medium width, natural meta pad</i>	21/22 - 48
971330-000		Multiforte, hautfarben, Shore ca. 40°, 13 mm <i>Multiforte, skin-colour, Shore approx. 40°, 13 mm</i>		lang, breit, Pelottendruck <i>full length, wide, natural meta pad</i>	19/20 - 50
977016-000		Thermofix, weiss, Shore ca. 45°, 14 mm <i>Thermofix, white, Shore approx. 45°, 14 mm</i>		Viertelschale, lang, Pelottendruck <i>quarter cup, full length, natural meta pad</i>	19/20 - 48
977021-000		Thermofix, braun, Shore ca. 45°, 14 mm <i>Thermofix, brown, Shore approx. 45°, 14 mm</i>		Viertelschale, lang, Pelottendruck <i>quarter cup, full length, natural meta pad</i>	19/20 - 48
977214-000		Thermofix, braun, Shore ca. 45°, 14 mm <i>Thermofix, brown, Shore approx. 45°, 14 mm</i>		Viertelschale, kurz, Pelottendruck <i>quarter cup, short, natural meta pad</i>	35 - 47
977016-006		Thermofix, weiss, Shore ca. 45°, 14 mm <i>Thermofix, white, Shore approx. 45°, 14 mm</i>		Viertelschale, lang, Pelottendruck, G-schliff <i>quarter cup, full length, natural meta pad, grinded bottom</i>	19/20 - 48
977016-300		Thermofix, weiss, Shore ca. 45°, 14 mm <i>Thermofix, white, Shore approx. 45°, 14 mm</i>		Viertelschale, lang, ohne Pelottendruck <i>quarter cup, full length, without meta pad</i>	19/20 - 48



470870-000



971330-000



977016-000



977021-000



977214-000

Mehrschicht
multi layer


470622-000



470754-000



977313-000



977710-000

Artikel-Nr. Code-No.	Unterbau/Kern Base	Decke Top Layer	Bezug Top Cover	Ausführung Design	Größen Sizes
470622-000	Multiform, weiss, Shore ca. 30°, 7 mm, LG-Keil <i>Multiform, white, Shore approx. 30°, 7 mm, longitudinal arch wedge</i>	Multiform, blau, Shore ca. 30°, 4 mm, perforiert <i>Multiform, blue, Shore approx. 30°, 4 mm, perforated</i>		lang, mittel, Pelottendruck <i>full length, medium width, natural meta pad</i>	19/20 - 48
470754-000	Multiform, weiss, Shore ca. 30°, 7 mm, LG-Keil <i>Multiform, white, Shore approx. 30°, 7 mm, longitudinal arch wedge</i>	Multiform, weiss, Shore ca. 30°, 4 mm, perforiert <i>Multiform, white, Shore approx. 30°, 4 mm, perforated</i>		lang, mittel, Pelottendruck <i>full length, medium width, natural meta pad</i>	19/20 - 48
977313-000	Multiform, weiss, Shore ca. 30°, 10 mm <i>Multiform, white, Shore approx. 30°, 10 mm</i>	Multiform, blau, Shore ca. 30°, 4 mm, perforiert <i>Multiform, blue, Shore approx. 30°, 4 mm, perforated</i>		Viertelschale, lang, Pelottendruck <i>quarter cup, full length, natural meta pad</i>	19/20 - 48
977710-000	Multiform, weiss, Shore ca. 30°, 10 mm <i>Multiform, white, Shore approx. 30°, 10 mm</i>	Multiform, blau, Shore ca. 30°, 4 mm, perforiert <i>Multiform, blue, Shore approx. 30°, 4 mm, perforated</i>		lang, breit, Pelottendruck <i>full length, wide, natural meta pad</i>	19/20 - 48
977710-095	Multiform, weiss, Shore ca. 30°, 10 mm <i>Multiform, white, Shore approx. 30°, 10 mm</i>	Multiform, hautfarben, Shore ca. 30°, 4 mm, perforiert <i>Multiform, skin-colour, Shore approx. 30°, 4 mm, perforated</i>		lang, breit, Pelottendruck <i>full length, wide, natural meta pad</i>	19/20 - 48
977710-244	Multiform, weiss, Shore ca. 30°, 10 mm <i>Multiform, white, Shore approx. 30°, 10 mm</i>	Multiform, weiss, Shore ca. 30°, 3 mm <i>Multiform, white, Shore approx. 30°, 3 mm</i>		lang, breit, Pelottendruck <i>full length, wide, natural meta pad</i>	19/20 - 48

Mehrschicht
multi layer

Artikel-Nr. Code-No.	Unterbau/Kern Base	Decke Top Layer	Bezug Top Cover	Ausführung Design	Größen Sizes
977735-000	Multiform, weiss, Shore ca. 30°, 10 mm, THK - Verstärkung <i>Multiform, white, Shore approx. 30°, 10 mm, THK reinforcement</i>	Multiform, blau, Shore ca. 30°, 4 mm, perforiert <i>Multiform, blue, Shore approx. 30°, 4 mm, perforated</i>		lang, mittel, Gummipolotte, vorne offen <i>full length, medium width, rubber pad, front open</i>	35 - 48
975010-000	Multiform, weiss, Shore ca. 30°, 3 mm, THK - Verstärkung <i>Multiform, white, Shore approx. 30°, 3 mm, THK reinforcement</i>	Multiform, weiss, Shore ca. 30°, 3 mm <i>Multiform, white, Shore approx. 30°, 3 mm</i>		lang, breit, Pelottendruck <i>full length, wide, natural meta pad</i>	35 - 48
977715-000	Multiform, weiss, Shore ca. 30°, 10 mm, THK - Verstärkung <i>Multiform, white, Shore approx. 30°, 10 mm, THK reinforcement</i>	Multiform, blau, Shore ca. 30°, 4 mm, perforiert <i>Multiform, blue, Shore ca. 30°, 4 mm, perforated</i>		lang, breit, Gummipolotte, vorne offen <i>full length, wide, rubber pad, front open</i>	19/20 - 48
977753-000	Multiform, weiss, Shore ca. 30°, 10 mm, THK - Verstärkung <i>Multiform, white, Shore approx. 30°, 10 mm, THK reinforcement</i>	Multiform, weiss, Shore ca. 30°, 7 mm <i>Multiform, white, Shore approx. 30°, 7 mm</i>		Naturform, lang, Pelottendruck <i>natural shape, full length, natural meta pad</i>	35 - 48
977838-000	Multiform, weiss, Shore ca. 30°, 10 mm, THK - Verstärkung <i>Multiform, white, Shore approx. 30°, 10 mm, THK reinforcement</i>	Multiform, weiss, Shore ca. 30°, 4 mm <i>Multiform, white, Shore approx. 30°, 4 mm</i>		lang, breit, Gummipolotte <i>full length, wide, rubber pad</i>	35 - 50
977745-164	Multiform, weiss, Shore ca. 30°, 10 mm, THK - Verstärkung <i>Multiform, white, Shore approx. 30°, 10 mm, THK reinforcement</i>	Multiform, braun, Shore ca. 30°, 4 mm, perforiert <i>Multiform, brown, Shore approx. 30°, 4 mm, perforated</i>		lang, mittel, Gummipolotte, vorne offen <i>full length, medium width, rubber pad, front open</i>	35 - 48
975010-300	Multiform, weiss, Shore ca. 30°, 3 mm, THK - Verstärkung <i>Multiform, white, Shore approx. 30°, 3 mm, THK reinforcement</i>	Multiform, weiss, Shore ca. 30°, 3 mm <i>Multiform, white, Shore approx. 30°, 3 mm</i>		lang, breit, ohne Pelottendruck <i>full length, wide, without meta pad</i>	35 - 48



975010-000



977715-000



977753-000



977838-000



977745-164

Mehrschicht
multi layer


470614-000



977841-000



976841-050



977713-000

Artikel-Nr. Code-No.	Unterbau/Kern Base	Decke Top Layer	Bezug Top Cover	Ausführung Design	Größen Sizes
470614-000	Multiform, weiss, Shore ca. 30°, 10 mm, THK - Versteifung <i>Multiform, white, Shore approx. 30°, 10 mm, THK reinforcement</i>	Multiform, blau, Shore ca. 30°, 5 mm, perforiert <i>Multiform, blue, Shore approx. 30°, 5 mm, perforated</i>		Schalenform, lang, Pelottendruck <i>cup, full length, natural meta pad</i>	19/20 - 48
976841-000	Multiform, weiss, Shore ca. 30°, 10 mm, THK - Versteifung <i>Multiform, white, Shore approx. 30°, 10 mm, THK reinforcement</i>	p ² , hautfarben, Shore ca. 18°, 3 mm <i>p², skin-colour, Shore approx. 18°, 3 mm</i>		vorne offen, lang, mittel, Gummipelotte <i>front open, full length, medium width, rubber pad</i>	19/20 - 48
977841-000	Multiform, weiss, Shore ca. 30°, 10 mm, THK - Versteifung <i>Multiform, white, Shore approx. 30°, 10 mm, THK reinforcement</i>	p ² , hautfarben, Shore ca. 18°, 3 mm <i>p², skin-colour, Shore approx. 18°, 3 mm</i>	Rindleder, perforiert, mild <i>neat's leather, perforated, soft</i>	Naturform, lang, Gummipelotte <i>natural shape, full length, rubber pad</i>	35 - 48
976841-035	Multiform, weiss, Shore ca. 30°, 10 mm, THK - Versteifung <i>Multiform, white, Shore approx. 30°, 10 mm, THK reinforcement</i>	p ² , hautfarben, Shore ca. 18°, 3 mm <i>p², skin-colour, Shore approx. 18°, 3 mm</i>	Rindleder, perforiert, mild <i>neat's leather, perforated, soft</i>	lang, mittel Gummipelotte <i>full length, medium width, rubber pad</i>	19/20 - 48
976841-050	Multiform, weiss, Shore ca. 30°, 10 mm, THK - Versteifung <i>Multiform, white, Shore approx. 30°, 10 mm, THK reinforcement</i>	p ² , hautfarben, Shore ca. 18°, 3 mm <i>p², skin-colour, Shore approx. 18°, 3 mm</i>	Rindleder, perforiert, mild <i>neat's leather, perforated, soft</i>	lang, mittel, Gummipelotte, vorne offen <i>full length, medium width, rubber pad, front open</i>	19/20 - 48
977713-156	Multiform, weiss, Shore ca. 30°, 10 mm, THK - Versteifung <i>Multiform, white, Shore approx. 30°, 10 mm, THK reinforcement</i>	Multiform, braun, Shore ca. 30°, 4 mm, perforiert <i>Multiform, brown, Shore approx. 30°, 4 mm, perforated</i>		Viertelschale, lang, Pelottendruck <i>quarter cup, full length, natural meta pad</i>	19/20 - 48
977713-243	Multiform, weiss, Shore ca. 30°, 10 mm, THK - Versteifung <i>Multiform, white, Shore approx. 30°, 10 mm, THK reinforcement</i>	Multiform, weiss, Shore ca. 30°, 4 mm <i>Multiform, white, Shore approx. 30°, 4 mm</i>		Viertelschale, lang, Pelottendruck <i>quarter cup, full length, natural meta pad</i>	19/20 - 48
977713-000	Multiform, weiss, Shore ca. 30°, 10 mm, THK - Versteifung <i>Multiform, white, Shore approx. 30°, 10 mm, THK reinforcement</i>	Multiform, blau, Shore ca. 30°, 4 mm, perforiert <i>Multiform, blue, Shore approx. 30°, 4 mm, perforated</i>		Viertelschale, lang, Pelottendruck <i>quarter cup, full length, natural meta pad</i>	19/20 - 48

Mehrschicht
multi layer

Artikel-Nr. Code-No.	Unterbau/Kern Base	Decke Top Layer	Bezug Top Cover	Ausführung Design	Größen Sizes
470617-000	Multiform, weiss, Shore ca. 30°, 7 mm, LG-Keil, THK - Verstärkung Multiform, white, Shore approx. 30°, 7 mm, longitudinal arch wedge, THK reinforcement	Multiform, braun, Shore ca. 30°, 7 mm, perforiert Multiform, brown, Shore approx. 30°, 7 mm, perforated		Naturform, lang, Pelottendruck natural shape, full length, natural meta pad	19/20 - 48
470624-000	Multiform, weiss, Shore ca. 30°, 7 mm, LG-Keil, THK - Verstärkung Multiform, white, Shore approx. 30°, 7 mm, longitudinal arch wedge, THK reinforcement	Multiform, blau, Shore ca. 30°, 4 mm, perforiert Multiform, blue, Shore approx. 30°, 4 mm, perforated		lang, mittel, Pelottendruck full length, medium width, natural meta pad	19/20 - 48
470778-000	Multiform, weiss, Shore ca. 30°, 7 mm, LG-Keil, THK - Verstärkung Multiform, white, Shore approx. 30°, 7 mm, longitudinal arch wedge, THK reinforcement	Multiform, weiss, Shore ca. 30°, 4 mm, perforiert Multiform, white, Shore approx. 30°, 4 mm, perforated		lang, mittel, Pelottendruck full length, medium width, natural meta pad	19/20 - 48
470810-000	Multiform, weiss, Shore ca. 30°, 7 mm, LG-Keil, THK - Verstärkung Multiform, white, Shore approx. 30°, 7 mm, longitudinal arch wedge, THK reinforcement	Multiform, multicolor, Shore ca. 30°, 4 mm, perforiert Multiform, multicolor, Shore approx. 30°, 4 mm, perforated		lang, mittel, Pelottendruck full length, medium width, natural meta pad	19/20 - 48
470723-000	Multiform, weiss, Shore ca. 30°, 7 mm, LG-Keil, THK - Verstärkung Multiform, white, Shore approx. 30°, 7 mm, longitudinal arch wedge, THK reinforcement	Multiform, hautfarben, Shore ca. 30°, 4 mm, perforiert Multiform, skin-colour, Shore approx. 30°, 4 mm, perforated		lang, mittel, Pelottendruck full length, medium width, natural meta pad	19/20 - 48
471730-000	Multiform, hautfarben, Shore ca. 30°, 7 mm, LG-Keil, THK - Verstärkung Multiform, skin-colour, Shore approx. 30°, 7 mm, longitudinal arch wedge, THK reinforcement	Multiform, braun, Shore ca. 30°, 4 mm, perforiert Multiform, brown, Shore approx. 30°, 4 mm, perforated		lang, mittel, Pelottendruck full length, medium width, natural meta pad	19/20 - 48
471732-000	Multiform, weiss, Shore ca. 30°, 7 mm, LG-Keil, THK - Verstärkung Multiform, white, Shore approx. 30°, 7 mm, longitudinal arch wedge, THK reinforcement	Multiform, hautfarben, Shore ca. 30°, 4 mm, perforiert Multiform, skin-colour, Shore approx. 30°, 4 mm, perforated		Naturform, lang, Pelottendruck natural shape, full length, natural meta pad	19/20 - 48




470617-000



470624-000



470778-000

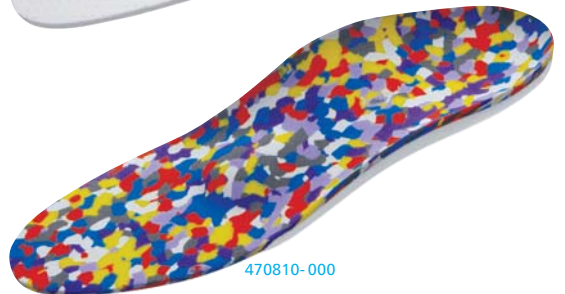

470810-000



470723-000



471732-000

Mehrschicht
multi layer


470730-000



470733-000



471625-000



471723-000

Artikel-Nr. Code-No.	Unterbau/Kern Base	Decke Top Layer	Bezug Top Cover	Ausführung Design	Größen Sizes
470730-000	Multiform, hautfarben, Shore ca. 30°, 7 mm, LG-Keil, THK - Versteifung <i>Multiform, skin-colour, Shore approx. 30°, 7 mm, longitudinal arch wedge, THK reinforcement</i>	Multiform, braun, Shore ca. 30°, 4 mm, perforiert <i>Multiform, brown, Shore approx. 30°, 4 mm, perforated</i>		lang, mittel, Gummipelotte <i>full length, medium width, rubber pad</i>	19/20 - 48
470733-000	Multiform, weiss, Shore ca. 30°, 7 mm, LG-Keil, THK - Versteifung <i>Multiform, white, Shore approx. 30°, 7 mm, longitudinal arch wedge, THK reinforcement</i>	Multiform, hautfarben, Shore ca. 30°, 4 mm, perforiert <i>Multiform, skin-colour, Shore approx. 30°, 4 mm, perforated</i>		Naturform, lang, Gummipelotte <i>natural shape, full length, rubber pad</i>	19/20 - 48
471608-000	Multiform, weiss, Shore ca. 30°, 7 mm, LG-Keil, THK - Versteifung <i>Multiform, white, Shore approx. 30°, 7 mm, longitudinal arch wedge, THK reinforcement</i>	Multiform, blau, Shore ca. 30°, 4 mm, perforiert <i>Multiform, blue, Shore approx. 30°, 4 mm, perforated</i>		lang, breit, Gummipelotte <i>full length, wide, rubber pad</i>	19/20 - 48
471613-000	Multiform, hautfarben, Shore ca. 30°, 7 mm, LG-Keil, THK - Versteifung <i>Multiform, skin-colour, Shore approx. 30°, 7 mm, longitudinal arch wedge, THK reinforcement</i>	Multiform, braun, Shore ca. 30°, 4 mm, perforiert <i>Multiform, brown, Shore approx. 30°, 4 mm, perforated</i>		lang, breit, Gummipelotte <i>full length, wide, rubber pad</i>	19/20 - 48
471625-000	Multiform, weiss, Shore ca. 30°, 7 mm, LG-Keil, THK - Versteifung <i>Multiform, white, Shore approx. 30°, 7 mm, longitudinal arch wedge, THK reinforcement</i>	Multiform, blau, Shore ca. 30°, 4 mm, perforiert <i>Multiform, blue, Shore approx. 30°, 4 mm, perforated</i>		lang, mittel, Gummipelotte <i>full length, medium width, rubber pad</i>	19/20 - 48
471723-000	Multiform, weiss, Shore ca. 30°, 7 mm, LG-Keil, THK - Versteifung <i>Multiform, white, Shore approx. 30°, 7 mm, longitudinal arch wedge, THK reinforcement</i>	Multiform, hautfarben, Shore ca. 30°, 4 mm, perforiert <i>Multiform, skin-colour, Shore approx. 30°, 4 mm, perforated</i>		lang, mittel, Gummipelotte <i>full length, medium width, rubber pad</i>	19/20 - 48
471733-000	Multiform, braun, Shore ca. 30°, 7 mm, LG-Keil, THK - Versteifung <i>Multiform, brown, Shore approx. 30°, 7 mm, longitudinal arch wedge, THK reinforcement</i>	Multiform, braun, Shore ca. 30°, 4 mm, perforiert <i>Multiform, brown, Shore approx. 30°, 4 mm, perforated</i>		Naturform, lang, Gummipelotte, vorne offen <i>natural shape, full length, rubber pad, front open</i>	19/20 - 48
471774-000	Multiform, weiss, Shore ca. 30°, 7 mm, LG-Keil, THK - Versteifung <i>Multiform, white, Shore approx. 30°, 7 mm, longitudinal arch wedge, THK reinforcement</i>	Multiform, weiss, Shore ca. 30°, 4 mm, perforiert <i>Multiform, white, Shore approx. 30°, 4 mm, perforated</i>		lang, breit, Gummipelotte <i>full length, wide, rubber pad</i>	19/20 - 48

Mehrschicht
multi layer

Artikel-Nr. Code-No.	Unterbau/Kern Base	Decke Top Layer	Bezug Top Cover	Ausführung Design	Größen Sizes
471776-000	Multiform, weiss, Shore ca. 30°, 7 mm, LG-Keil, THK - Verstärkung <i>Multiform, white, Shore approx. 30°, 7 mm, longitudinal arch wedge, THK reinforcement</i>	Multiform, weiss, Shore ca. 30°, 4 mm, perforiert <i>Multiform, white, Shore approx. 30°, 4 mm, perforated</i>		lang, mittel, Gummipelotte <i>full length, medium width, rubber pad</i>	19/20 - 48
471778-000	Multiform, weiss, Shore ca. 30°, 7 mm, LG-Keil, THK - Verstärkung <i>Multiform, white, Shore approx. 30°, 7 mm, longitudinal arch wedge, THK reinforcement</i>	Multiform, weiss, Shore ca. 30°, 4 mm, perforiert <i>Multiform, white, Shore approx. 30°, 4 mm perforated</i>		lang, mittel, Gummipelotte <i>full length, medium width, rubber pad</i>	19/20 - 48
471778-104	Multiform, weiss, Shore ca. 30°, 7 mm, LG-Keil, THK - Verstärkung <i>Multiform, white, Shore approx. 30°, 7 mm, longitudinal arch wedge, THK reinforcement</i>	Multiform, weiss, Shore ca. 30°, 4 mm, perforiert <i>Multiform, white, Shore approx. 30°, 4 mm, perforated</i>	Ziegennarben- spalt, perforiert <i>goat leather, perforated</i>	lang, mittel, Gummipelotte <i>full length, medium width, rubber pad</i>	19/20 - 48
470709-000	Multiform, hautfarben, Shore ca. 30°, 7 mm, LG-Keil, THK - Verstärkung <i>Multiform, skin-colour, Shore approx. 30°, 7 mm, longitudinal arch wedge, THK reinforcement</i>	Multiform, braun, Shore ca. 30°, 4 mm, perforiert <i>Multiform, brown, Shore approx. 30°, 4 mm, perforated</i>		lang, mittel, retrokapital, Pelottendruck <i>full length, medium width, retrocapital support, natural meta pad</i>	35 - 48
977858-000	PE-Schaum, grau, Shore ca. 20°, 5 mm, LG-Keil, THK - Verstärkung <i>PE-Foam, grey, Shore approx. 20°, 5 mm, longitudinal arch wedge, THK reinforcement</i>	PE-Schaum, grau, Shore ca. 20°, 5 mm <i>PE-Foam, grey, Shore approx. 20°, 5 mm</i>		Naturform, lang, ohne Pelottendruck <i>natural shape, full length, without meta pad</i>	35 - 48
470791-000	PE-Schaum, grau, Shore ca. 20°, 5 mm, LG-Keil, THK - Verstärkung <i>PE-Foam, grey, Shore approx. 20°, 5 mm, longitudinal arch wedge, THK reinforcement</i>	PE-Schaum, grau, Shore ca. 20°, 5 mm <i>PE-Foam, grey, Shore approx. 20°, 5 mm</i>	Rindleder, perforiert, mild <i>neat's leather, perforated, soft</i>	lang, mittel, Gummipelotte <i>full length, medium width, rubber pad</i>	35 - 48
977873-133	PE-Schaum, grau, Shore ca. 20°, 5 mm, LG-Keil, THK - Verstärkung <i>PE-Foam, grey, Shore approx. 20°, 5 mm, longitudinal arch wedge, THK reinforcement</i>	PE-Schaum, grau, Shore ca. 20°, 5 mm <i>PE-Foam, grey, Shore approx. 20°, 5 mm</i>	Rindleder, perforiert, mild <i>neat's leather, perforated, soft</i>	lang, mittel, Gummipelotte <i>full length, medium width, rubber pad</i>	35 - 48



471776-000



471778-000



470709-000



977858-000



470791-000

NovaPED soft

Mehrschicht
multi layer



470610-000



471885-000



474097-000



977426-000



474156-000



471320-057

Artikel-Nr. Code-No.	Unterbau/Kern Base	Decke Top Layer	Bezug Top Cover	Ausführung Design	Größen Sizes
470610-000	Multiforte, weiss, Shore ca. 40°, 10 mm, THK - Verstärkung <i>Multiforte, white, Shore approx. 40°, 10 mm, THK reinforcement</i>	Multiform, weiss, Shore ca. 30°, 4 mm, perforiert <i>Multiform, white, Shore approx. 30°, 4 mm, perforated</i>		lang, mittel, Pelottendruck <i>full length, medium width, natural meta pad</i>	19/20 - 48
471885-000	Multiforte, blau, Shore ca. 40°, 12 mm <i>Multiforte, blue, Shore approx. 40°, 12 mm</i>	Multiform, weiss, Shore ca. 30°, 4 mm, perforiert <i>Multiform, white, Shore approx. 30°, 4 mm, perforated</i>		lang, breit, Pelottendruck, Unterseite flach <i>full length, medium width, natural meta pad, flat bottom</i>	35 - 50
474097-000	Multiforte, blau, Shore ca. 40°, 12 mm, LG-Keil, THK - Verstärkung <i>Multiforte, blue, Shore approx. 40°, 12 mm, longitudinal arch wedge, THK reinforcement</i>	Multicom, blau/blau, Shore ca. 30°/18°, 4 mm, perforiert <i>Multicom, blue/blue, Shore approx. 30°/18°, 4 mm, perforated</i>		Naturform, lang, Pelottendruck <i>natural shape, full length, natural meta pad</i>	35 - 46
977426-000	Multiforte, weiss, Shore ca. 40°, 10 mm, THK - Verstärkung <i>Multiforte, white, Shore approx. 40°, 10 mm, THK reinforcement</i>	Multiform, weiss, Shore ca. 30°, 5 mm <i>Multiform, white, Shore approx. 30°, 5 mm</i>		lang, breit, Pelottendruck <i>full length, wide, natural meta pad</i>	19/20 - 48
474156-000	Multiforte, blau, Shore ca. 40°, 7 mm, LG-Keil <i>Multiforte, blue, Shore approx. 40°, 7 mm, longitudinal arch wedge</i>	Multiform, blau, Shore ca. 30°, 4 mm <i>Multiform, blue, Shore approx. 30°, 4 mm</i>		lang, mittel, Pelottendruck <i>full length, medium width, natural meta pad</i>	35 - 48
470883-020	Multiforte, blau, Shore ca. 40°, 12 mm <i>Multiforte, blue, Shore approx. 40°, 12 mm</i>	Multicom, blau/blau, Shore ca. 30°/18°, 4 mm, perforiert <i>Multicom, blue/blue, Shore approx. 30°/18°, 4 mm, perforated</i>		lang, breit, Pelottendruck, Unterseite flach <i>full length, wide, natural meta pad, flat bottom</i>	35 - 48
471320-057	Multiforte, weiss, Shore ca. 40°, 10 mm <i>Multiforte, white, Shore approx. 40°, 10 mm</i>	ARUfoam 20, weiss, Shore ca. 35°, 5 mm <i>ARUfoam 20, white, Shore approx. 35°, 5 mm</i>		lang, mittel, Pelottendruck <i>full length, medium width, natural meta pad</i>	19/20 - 48
471320-102	Multiforte, weiss, Shore ca. 40°, 10 mm <i>Multiforte, white, Shore approx. 40°, 10 mm</i>	Multiform, weiss, Shore ca. 30°, 4 mm <i>Multiform, white, Shore approx. 30°, 4 mm</i>		lang, mittel, Pelottendruck <i>full length, medium width, natural meta pad</i>	19/20 - 48

NovaPED soft

Mehrschicht
multi layer



Artikel-Nr. Code-No.	Unterbau/Kern Base	Decke Top Layer	Bezug Top Cover	Ausführung Design	Größen Sizes
470610- 156	Multiforte, weiss, Shore ca. 40°, 10 mm, THK - Versteifung	Multiform, braun, Shore ca. 30°, 4 mm, perforiert		lang, mittel, Pelottendruck	19/20 - 48
	Multiforte, white, Shore approx. 40°, 10 mm, THK reinforcement	Multiform, brown, Shore approx. 30°, 4 mm, perforated		full length, medium width, natural meta pad	
977364- 209	Multiforte, hautfarben, Shore ca. 40°, 13 mm	Multiform, blau, Shore ca. 30°, 4 mm, perforiert		lang, mittel, Gummipelotte	19/20 - 48
	Multiforte, skin-colour, Shore approx. 40°, 13 mm	Multiform, blue, Shore approx. 30°, 4 mm, perforated		full length, medium width, rubber pad	
471320- 223	Multiforte, weiss, Shore ca. 40°, 10 mm	Multiform, weiss, Shore ca. 30°, 4 mm,	Kalbfutterleder, perforiert	lang, mittel, Pelottendruck	19/20 - 48
	Multiforte, white, Shore approx. 40°, 10 mm	Multiform, white, Shore approx. 30°, 4 mm	calfskin, perforated	full length, medium width, natural meta pad	
471320- 468	Multiforte, weiss, Shore ca. 40°, 10 mm	Multiform, weiss, Shore ca. 30°, 4 mm	Alcantara, beige	lang, mittel, Pelottendruck	19/20 - 48
	Multiforte, white, Shore approx. 40°, 10 mm	Multiform, white, Shore approx. 30°, 4 mm	Alcantara, beige	full length, medium width, natural meta pad	
471626- 477	Multiforte, blau, Shore ca. 40°, 7 mm, LG-Keil, THK - Versteifung	Multicom, blau/blau Shore ca. 30°/18°, 4 mm, perforiert	Alcantara, schwarz	lang, mittel, Gummipelotte	19/20 - 48
	Multiforte, blue, Shore approx. 40°, 7 mm, longitudinal arch wedge, THK reinforcement	Multicom, blue/blue, Shore approx. 30°/18°, 4 mm, perforated	Alcantara, black	full length, medium width, rubber pad	

NovaPED soft

Mehrschicht
multi layer



972778-000

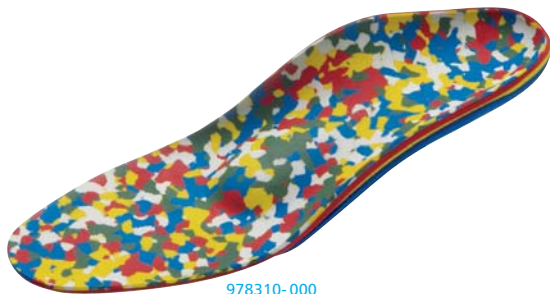

978310-000



470802-000



470807-000



972650-104

Artikel-Nr. Code-No.	Unterbau/Kern Base	Decke Top Layer	Bezug Top Cover	Ausführung Design	Größen Sizes
974325-000	Thermofix, braun, Shore ca. 45°, 10 mm <i>Thermofix, brown, Shore approx. 45°, 10 mm</i>	Pryx 25, braun, Shore ca. 25°, 6 mm <i>Pryx 25, brown, Shore approx. 25°, 6 mm</i>		Naturform, Viertelschale, lang, Pelottendruck <i>natural shape, quarter cup, full length, natural meta pad</i>	35 - 48
975013-000	Thermofix, weiss, Shore ca. 45°, 10 mm <i>Thermofix, white, Shore approx. 45°, 10 mm</i>	Multiform, weiss, Shore ca. 30°, 3 mm, perforiert <i>Multiform, white, Shore approx. 30°, 3 mm, perforated</i>		Viertelschale, lang, Pelottendruck <i>quarter cup, full length, natural meta pad</i>	35 - 48
972778-000	Pryx, antistatisch, schwarz, Shore ca. 55°, 10 mm <i>Pryx, anti-static, black, Shore approx. 55°, 10 mm</i>	Pryx, antistatisch, schwarz, Shore ca. 55°, 6 mm <i>Pryx, anti-static, black, Shore approx. 55°, 6 mm</i>		Naturform, Viertelschale, lang, Pelottendruck <i>natural shape, quarter cup, full length, natural meta pad</i>	35 - 48
978310-000	EVA, blau/gelb, Shore ca. 50°/50°, 3 mm/3 mm <i>EVA, blue/yellow, Shore approx. 50°/50°, 3 mm/3 mm</i>	EVA, rot, Shore ca. 50°, 3 mm <i>EVA, red, Shore approx. 50°, 3 mm</i>	EVA, multicolor, Shore ca. 30°, 2,5 mm <i>EVA, multicolor, Shore approx. 30°, 2,5 mm</i>	lang, mittel, Pelottendruck <i>full length, medium width, natural meta pad</i>	35 - 48
470802-000	Multikork, Shore ca. 60°, 10 mm, THK - Verstärkung <i>Multikork, blue, Shore approx. 60°, 10 mm, THK reinforcement</i>	Multiform, hautfarben, Shore ca. 30°, 4 mm, perforiert <i>Multiform, skin-colour, Shore approx. 30°, 4 mm, perforated</i>		lang, mittel, Pelottendruck <i>full length, medium width, natural meta pad</i>	19/20 - 48
470804-000	Multikork, Shore ca. 60°, 10 mm, THK - Verstärkung <i>Multikork, Shore approx. 60°, 10 mm, THK reinforcement</i>	Multiform, blau, Shore ca. 30°, 4 mm, perforiert <i>Multiform, blue, Shore approx. 30°, 4 mm, perforated</i>		lang, mittel, Pelottendruck <i>full length, medium width, natural meta pad</i>	19/20 - 48
470807-000	Multikork, Shore ca. 60°, 10 mm, THK - Verstärkung <i>Multikork, Shore approx. 60°, 10 mm, THK reinforcement</i>	Multiform, blau, Shore ca. 30°, 5 mm, perforiert <i>Multiform, blue, Shore approx. 30°, 5 mm, perforated</i>		Schalenform, lang, Pelottendruck <i>cup, full length, natural meta pad</i>	19/20 - 48
972650-104	Multiform, weiss, fest, Shore ca. 50°, 10 mm <i>Multiform, white, hard, Shore approx. 50°, 10 mm</i>	Multiform, weiss, Shore ca. 30°, 4 mm <i>Multiform, white, Shore approx. 30°, 4 mm</i>	Ziegennarben- spalt, perforiert <i>goat leather, perforated</i>	Naturform, lang, Pelottendruck <i>natural shape, full length, natural meta pad</i>	35 - 48

Kindergrößen bis 33/34 ohne Pelotten bzw. Pelottendruck. Children's sizes up to 33/34 without pad respectively natural meta pad.

Artikel-Nr. Code-No.	Unterbau/Kern Base	Decke Top Layer	Bezug Top Cover	Ausführung Design	Größen Sizes
972857-000	EVA, hautfarben, Shore ca. 25°, 6 mm, Spezial-Kern <i>EVA, skin-colour, Shore approx. 25°, 6 mm, special core</i>	EVA, hautfarben, Shore ca. 18°, 6 mm <i>EVA, skin-colour, Shore approx. 18°, 6 mm</i>	scheinova, beige <i>scheinova, beige</i>	Naturform, lang, Pelottendruck <i>natural shape, full length, natural meta pad</i>	35 - 48
972870-000	EVA, hautfarben, Shore ca. 25°, 6 mm, Spezial-Kern <i>EVA, skin-colour, Shore approx. 25°, 6 mm, special core</i>	EVA, hautfarben, Shore ca. 18°, 6 mm <i>EVA, skin-colour, Shore approx. 18°, 6 mm</i>		Naturform, lang, ohne Pelottendruck <i>natural shape, full length, without meta pad</i>	35 - 48
470899-000	Multiforte, blau, Shore ca. 40°, 5 mm, LG-Keil, THK - Verstärkung <i>Multiforte, blue, Shore approx. 40°, 5 mm, longitudinal arch wedge, THK reinforcement</i>	Multicom, blau/blau, Shore ca. 30°/18°, 4 mm, perforiert <i>Multicom, blue/blue, Shore approx. 30°/18°, 4 mm, perforated</i>		Naturform, lang, Pelottendruck <i>natural shape, full length, natural meta pad</i>	35 - 48
470662-000	Multiforte, weiss, Shore ca. 40°, 5 mm, LG-Keil, THK - Verstärkung <i>Multiforte, white, Shore approx. 40°, 5 mm, longitudinal arch wedge, THK reinforcement</i>	Plastocom, hautfarben/blau, Shore ca. 20°/18°, 5 mm <i>Plastocom, skin-colour/blue, Shore approx. 20°/18°, 5 mm</i>		Naturform, lang, Pelottendruck <i>natural shape, full length, natural meta pad</i>	19/20 - 48
470898-000	Multiforte, weiss, Shore ca. 40°, 5 mm, LG-Keil, THK - Verstärkung <i>Multiforte, white, Shore approx. 40°, 5 mm, longitudinal arch wedge, THK reinforcement</i>	Plastocom, hautfarben/blau, Shore ca. 20°/18°, 5 mm <i>Plastocom, skin-colour/blue, Shore approx. 20°/18°, 5 mm</i>		Naturform, lang, Pelottendruck <i>natural shape, full length, natural meta pad</i>	19/20 - 48
973545-000	Multiform, weiss, fest, Shore ca. 50°, 10 mm, THK - Verstärkung <i>Multiform, white,hard, Shore approx. 50°, 10 mm, THK reinforcement</i>	Multiform, blau, Shore ca. 30°, 4 mm <i>Multiform, blue, Shore approx. 30°, 4 mm</i>	ARUfoam, lachs, Shore ca. 14°, 3 mm <i>ARUfoam, salmon, Shore approx. 14°, 3 mm</i>	Naturform, lang, ohne Pelottendruck, vorne offen <i>natural shape, full length, without meta pad, front open</i>	35 - 48



972857-000



972870-000



470899-000



470662-000



470898-000



973545-000

